


## Capsular Tension Ring with Scleral Fixation Arm

This is a modified CTR which has an additional fixation eyelet that allows it to be sutured to the scleral wall without piercing the capsular bag. This ring is used in cases where a standard CTR may not provide adequate centration and stability.

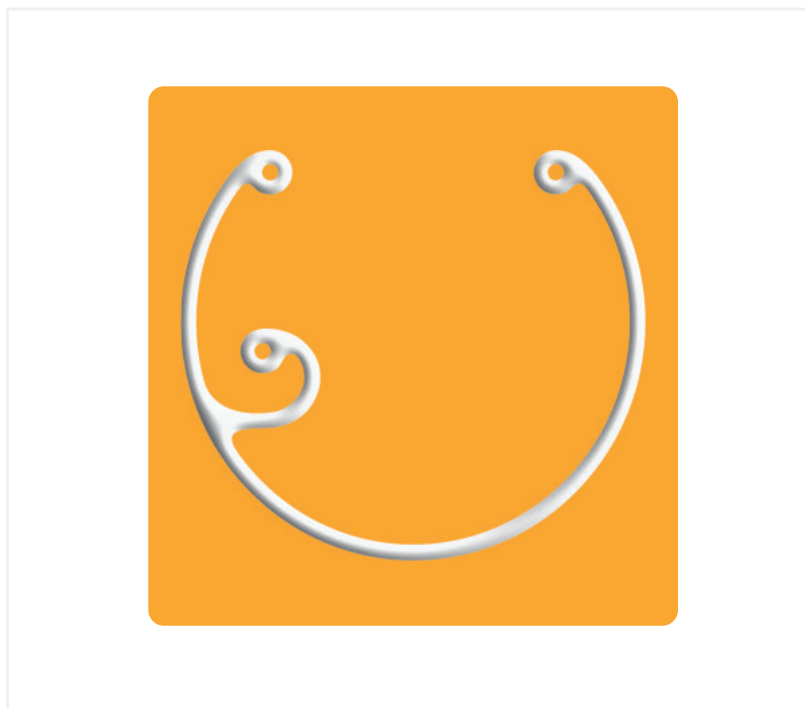

### Product Description

| Product Code | Model           | Design | Normal Size | Compressed Size | Eyelets | Material | Qty. per box    |
|--------------|-----------------|--------|-------------|-----------------|---------|----------|-----------------|
| MIPL/D3/T6   | CTR Type 6 (L)  |        | 12.0 mm     | 10.0 mm         | 3       | PMMA     | 1 pc. (Sterile) |
| MIPL/D3/T7   | CTR Type 7 (R)  |        | 12.0 mm     | 10.0 mm         | 3       | PMMA     | 1 pc. (Sterile) |
| MIPL/D3/T8   | CTR Type 8 (LR) |        | 12.0 mm     | 10.0 mm         | 4       | PMMA     | 1 pc. (Sterile) |

### Advantages

- It is used to stabilize and recenter the capsular bag
- Specially designed for scleral fixation with sutures through eyelets in a cataract surgery.
- The eyelet is projecting anteriorly which allows it to sit outside the capsular bag
- A double eyelet design is also available for cases which require additional support

## Indications

- Cases with extreme zonular weakness
- Cases with missing zonules larger than 4 clock hours
- Lens subluxation

## Instructions for Use

The instructions for use and surgical technique are supplied in the form of a leaflet with instructions given in diagram form where appropriate . The international symbols used on the packaging are also explained.

## Sterilization

Products are sterilized through Ethylene Oxide (EO) according to a validated cycle.



**MADHU INSTRUMENTS PVT LTD**

A-260, Okhla Industrial Area, Phase - 1, New Delhi - 110020, India

T: (+91-11) 4270 1030, 31 | E: [sales@madhuinstruments.com](mailto:sales@madhuinstruments.com) | [www.madhuinstruments.com](http://www.madhuinstruments.com)

Customer Care No. +91 9810866176 (for international queries), +91 7428186111 (for domestic queries)